



C-MAX

MOBILE ROBUST SAFE





THE C-MAX IN DETAIL

Total weight, approved	172 kg
Max. weight of the person	140 kg
Climbing speed ⁽¹⁾ (continually adjustable)	9 - 23 steps per minute ⁽¹⁾
Range with one charge of the batteries at an average person weight of 75 kg	approx. 25 levels ⁽¹⁾ (equates 375 steps)
Max. height of step	210 225 mm (extended)
Width	491 444 mm (without arm rests)
Depth	992 780 mm (folded in foot rest)
Height	1089 1431 mm
Batterys	2 x 12 V / 5 Ah
Direct current motor	24 V / 275 W
Weight, climbing unit	17.8 kg
Weight, C-Part with seat	5.3 kg
Weight, back rest incl. grip	4.5 kg
Weight, battery pack	4.0 kg
Weight, arm rests	1.3 kg
Total weight (incl. battery pack)	32.9 kg

(1) The climbing speed as well as the range with one charge of the batteries is reduced if the ComfortStep is activated. The climbing speed and the range also depend on



AAT Alber Antriebstechnik GmbH
Postfach 10 05 60 · D-72426 Albstadt

Phone: +49.7431.1295-0 · Fax: 1295-35
info@aat-online.de · www.aat-online.de



OVERCOMING STAIRS TROUGH

UNRESTRICTED MOBILITY

Your mobility is a precious commodity and significantly influences your quality of life. This is why it is essential to keep your freedom of movement as unrestricted as possible.

Our solution is called c-max - the mobile stair climbing aid that provides you with unlimited

mobility, not only at home, but also on trips!

No matter if you want to overcome stairs in your own home or explore new places - the c-max is always at your side. It enables you to master stairs effortlessly, no matter where you are.

The c-max is especially maneuverable and can be used without any difficulties on different stair surfaces like wood, stone, metal, PVC or carpet. The big wheels and robust climbing legs even allow to climb grating stairs.



**THE UNIVERSAL STAIR
CLIMBER**

WITH THE C-MAX

BENEFIT TO THE MAXIMUM



COMFORT STEP FUNCTION

For a comfortable ride and a soft touchdown of the climbing legs on every step. every step.



AUTOMATIC BRAKING FUNCTION

The automatic braking function when descending is part of our part of our sophisticated safety concept.




CAN BE USED AS A PUSH WHEELCHAIR

Thanks to its robust wheels, the c-max can also be used as a push wheelchair.




booklet_c-max_gb_12_2023 (1.1)

This brochure can also be displayed in big letters for easier readability. It is available for free download at www.aat-online.de. In addition, you will find product-specific user manuals with more technical data and further information on our homepage.

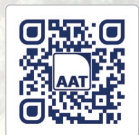

AAT online

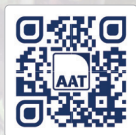

Contact agents



Social Media



The add-on kits of AAT coincide with the guidelines of the European Union 2006/42/EG, 2007/47/EG and 93/42/EG.

In the interest of progress we reserve the right for technical changes

AAT is certified according to DIN EN ISO 9001:2015 and DIN EN ISO 13485:2016.



Engineered in germany

